



- SS909HFC** Safety Station with WideArea™ Eye/Face Wash, Hand and Foot Control, Stainless Steel Bowl
- SS909P-HFC** Safety Station with WideArea™ Eye/Face Wash, Hand and Foot Control, Plastic Bowl

Application: Combination WideArea™ eye/face wash and shower safety station. Eye/face wash is activated by flag handle or foot treadle and features four GS-Plus™ spray-type outlet heads that deliver a flood of water over a wide area of coverage for complete rinsing of eyes and face.

Shower Head: 10" diameter orange ABS plastic with 20 GPM flow control.

Shower Valve: 1" IPS chrome plated brass stay-open ball valve. Valve is US-made with chrome plated brass ball and PTFE seals. Furnished with stainless steel actuating arm and 29" stainless steel pull rod.

Spray Head Assembly: Four GS-Plus™ spray heads. Each head has a "flip top" dust cover, internal flow control and filter to remove impurities from water flow.

Eye/Face Wash Bowl: 11-1/8" diameter stainless steel (SS909HFC) or 11-3/4" diameter orange ABS plastic (SS909P-HFC).

Eye/Face Wash Valve: 1/2" IPS chrome plated brass stay-open ball valve. Valve is US-made with chrome plated brass ball and PTFE seals. Valve is activated by flag handle or foot treadle. Unit remains in operation until handle is returned to closed position.

Pipe And Fittings: Schedule 40 galvanized steel. Furnished with orange polyethylene pipe covers for high visibility and corrosion resistance.

Supply: 1-1/4" NPT female top or side inlet.

Waste: 1-1/4" NPT female outlet.

Sign: ANSI-compliant identification sign.

Quality Assurance: Valve and spray head assemblies are factory assembled and water tested prior to shipment.

Available Options


- TM6040** Thermostatic mixing valve precisely blends hot and cold water to deliver tepid water as required by ANSI Z358.1-2014. Refer to "Thermostatic Mixing Valves" section for complete technical and product information.
- GC** Powder coated finish on galvanized pipe and fittings. Available colors include orange, yellow, red and green.
- BC** Stainless steel cover for eye/face wash bowl.
- GRN** Green ABS plastic shower head.
- YEL** Yellow ABS plastic shower head.
- SSH** Stainless steel shower head.
- AP275-200** Electric flashing light and alarm horn unit for mounting on wall or vertical pipe. Light illuminates and horn sounds when eye/face wash or shower is activated. Furnished complete with flow switch and mounting hardware.
- AP250-015** Modesty curtain for mounting on safety station.





□ **SS909HFC** Safety Station with WideArea™ Eye/Face Wash, Hand and Foot Control, Stainless Steel Bowl

□ **SS909P-HFC** Safety Station with WideArea™ Eye/Face Wash, Hand and Foot Control, Plastic Bowl

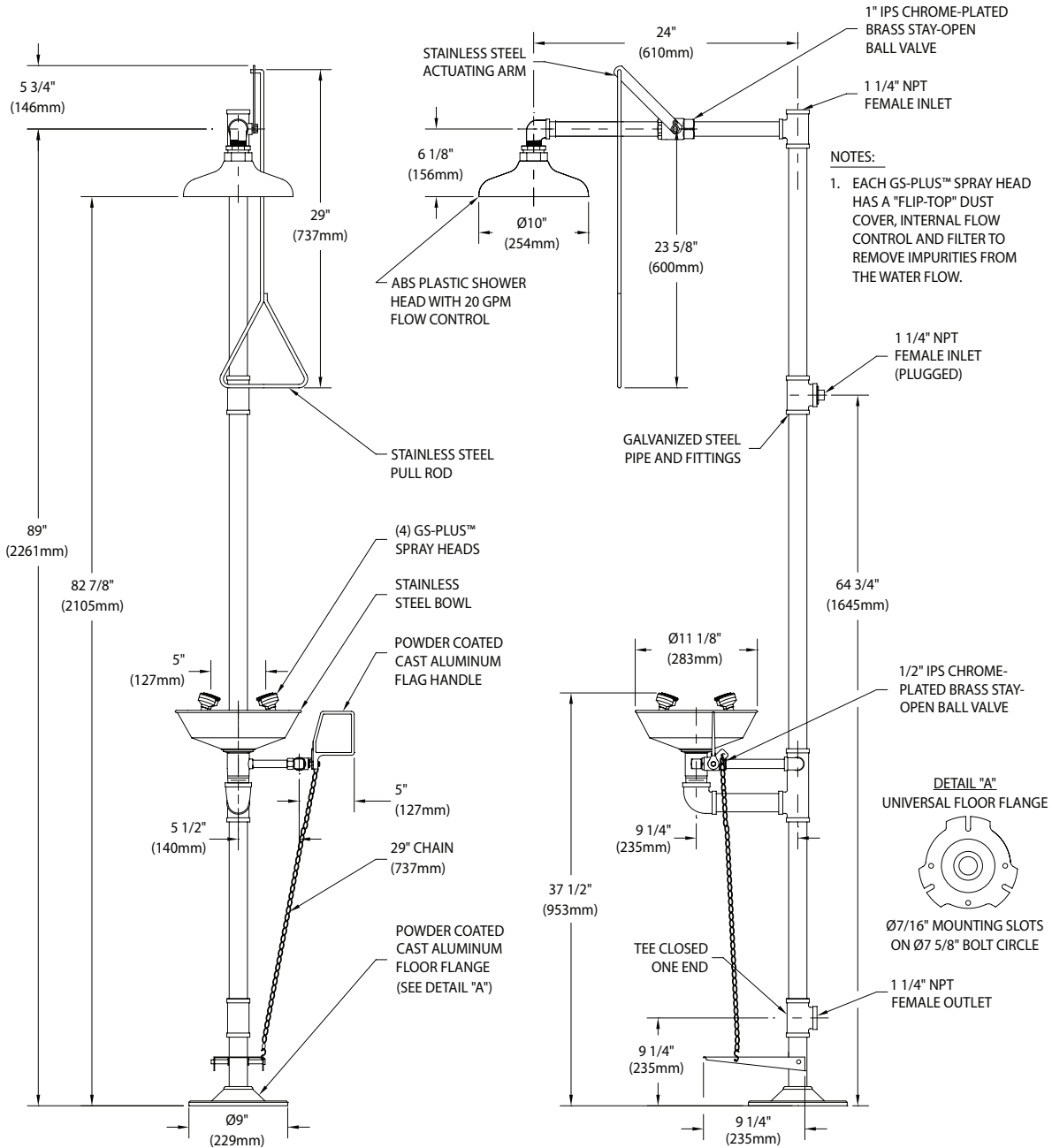

Rev. 022618

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Sign Included



Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are ± 1/4" (6mm). rev. 0522

WaterSaver Faucet Co. 312 666 5500 TELEPHONE
701 W Erie St 312 666 5501 FACSIMILE
Chicago, IL 60654 wsflab.com

Listed 8116. Units have been tested to and comply with ANSI Z358.1-2014 and the Uniform Plumbing Code.

