

BT1614BH-WSA Mixer Tap, Bench Mounted, Wrist Blade Handles

BT1614BHR-WSA Mixer Tap, Bench Mounted, Wrist Blade Handles (Reversed)



Application: Laboratory water mixer tap for hot and cold water.

Mounting: Bench or countertop mounted.

Body: Cast brass.

Headworks: WaterSaver self-contained compression valve cartridges with replaceable stainless steel seats.


Handles: Forged brass wrist blade handles with color coded index discs. Index discs are marked per EN 13792. On BT1614BH-WSA, handles rotate forward toward user to open, back to close. On BT1614BHR-WSA, handles rotate back to open, forward to close.

Gooseneck: 150mm (6") convertible rigid/swing gooseneck.


Inlets: (2) G1/2 male inlets on 180mm center dimension. Furnished with G1/2 mounting shanks, locknuts and washers.

Outlet: G3/8 female outlet with removable anti-splash serrated hose end.

Quality Assurance: Fitting is fully assembled and factory tested prior to shipment. 10 bar (150 PSI) test pressure and maximum working pressure.

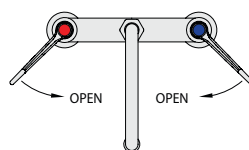

BT1614BH-WSA

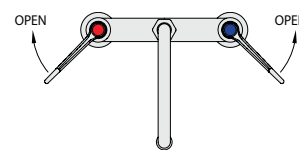

BT1614BHR-WSA

Ordering Information

Options		BT1614BH-WSA	BT1614BHR-WSA
Valve Construction	Compression valve cartridges	<input type="checkbox"/> BT1614BH-WSA	<input type="checkbox"/> BT1614BHR-WSA
	Ceramic disc valve cartridges	<input type="checkbox"/> BTC1614BH-WSA	<input type="checkbox"/> BTC1614BHR-WSA
Gooseneck	200mm (8") gooseneck in place of 150mm (6")	<input type="checkbox"/> BT1614-8BH-WSA	<input type="checkbox"/> BT1614-8BHR-WSA
Outlet	BO055 aerator in place of hose end	<input type="checkbox"/> BT1624WSA	<input type="checkbox"/> BT1624R-WSA
	BT118 pipe interrupter fitting with hose end	<input type="checkbox"/> BT1614-118BH-WSA	<input type="checkbox"/> BT1614-118BHR-WSA


Product Finish

<input type="checkbox"/> WHT	Gloss white powder coated finish	<input type="checkbox"/> PCL	Polished chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRY	Gloss light gray powder coated finish	<input type="checkbox"/> SCC	Satin chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRD	Gloss dark gray powder coated finish	<input type="checkbox"/> SNC	Satin nickel plated finish with clear epoxy coating
<input type="checkbox"/> STM	Starburst metallic powder coated finish		

Note: Refer to "Product Finishes" catalog for complete information on finishes.

Media

- BT1614BH-WSA** Mixer Tap, Bench Mounted, Wrist Blade Handles
- BT1614BHR-WSA** Mixer Tap, Bench Mounted, Wrist Blade Handles (Reversed)

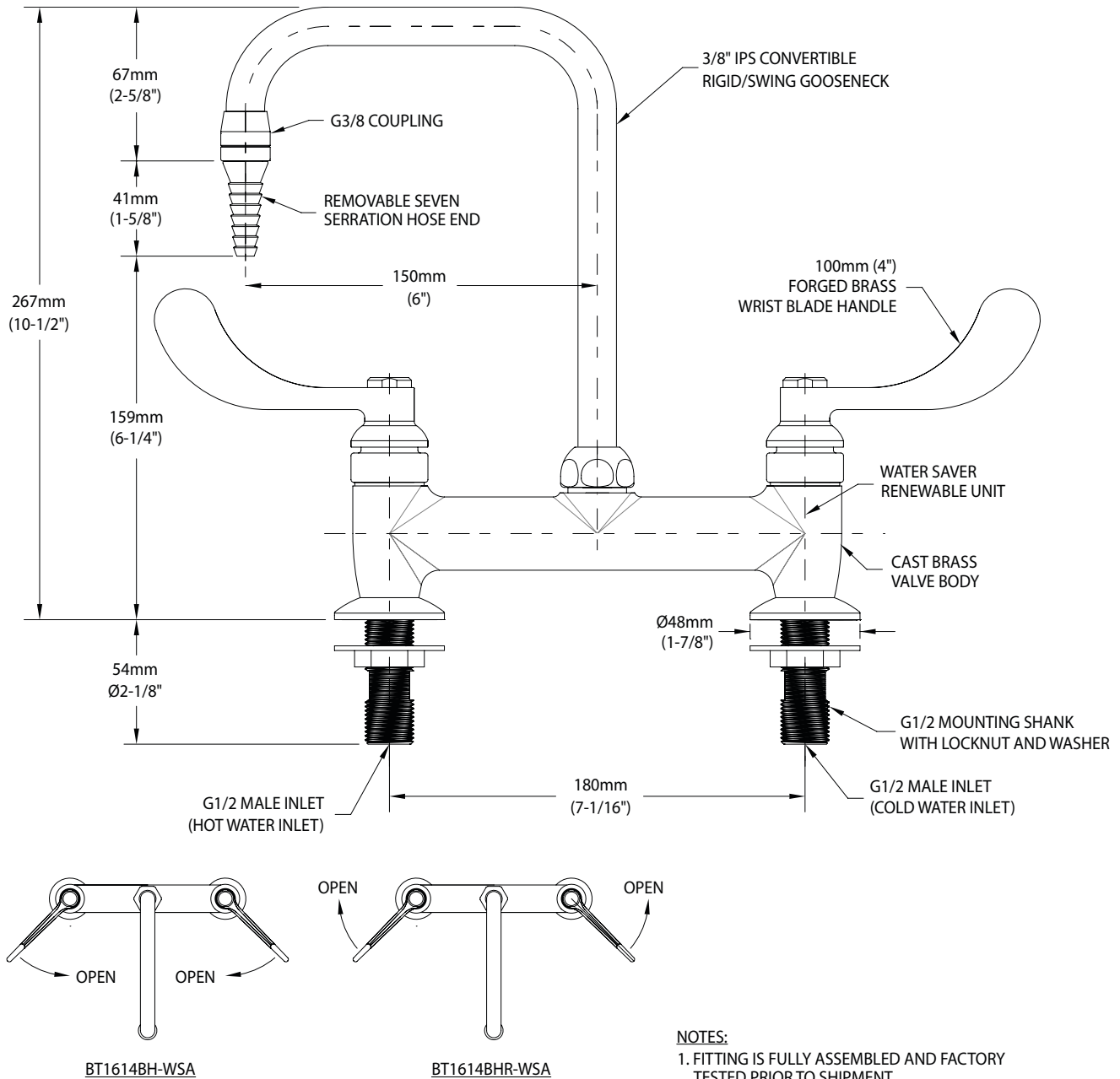

- NOTES:**
1. FITTING IS FULLY ASSEMBLED AND FACTORY TESTED PRIOR TO SHIPMENT.
 2. MAXIMUM WORKING PRESSURE IS 10 BAR (150 PSI).

Rev. 052318

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are ± 6mm (1/4"). Rev. 071718