

- SS902HS** Safety Station with Eyewash, Auxiliary Drench Hose, Stainless Steel Bowl
- SS902P-HS** Safety Station with Eyewash, Auxiliary Drench Hose, Plastic Bowl

**Application:** Combination eyewash and shower safety station. Furnished with auxiliary drench hose for rinsing user's eyes, face or body. Two GS-Plus™ spray-type outlet heads deliver a flood of water for rinsing eyes.

**Shower Head:** 10" diameter orange ABS plastic.

**Shower Valve:** 1" IPS chrome plated brass stay-open ball valve. Valve is US-made with chrome plated brass ball and Teflon® seals. Furnished with stainless steel actuating arm and 29" stainless steel pull rod.

**Spray Head Assembly:** Two GS-Plus™ spray heads. Each head has a "flip top" dust cover, internal flow control and filter to remove impurities from the water flow.

**Eyewash Bowl:** 11-1/2" diameter. Bowl is stainless steel (SS902HS) or orange ABS plastic (SS902P-HS).

**Eyewash Valve:** 1/2" IPS chrome plated brass stay-open ball valve. Valve is US-made with chrome plated brass ball and Teflon® seals.

**drench hose unit:** Auxiliary hand-held drench hose for rinsing eyes, face or body. Unit has single GS-Plus™ spray type outlet head, self-closing valve with squeeze handle, 8 ft. reinforced PVC hose and wall hook.

**Pipe and Fittings:** Schedule 40 galvanized steel. Furnished with orange polyethylene covers for vertical piping for high visibility and corrosion resistance.

**Supply:** 1-1/4" NPT female top or side inlet.

**Waste:** 1-1/4" NPT female outlet. Outlet can be positioned at either 9-1/4" or 19-5/8" above finished floor by reversing lower pipe nipples.

**Sign:** Furnished with ANSI-compliant identification sign.

**Quality Assurance:** Valve and spray head assemblies are factory assembled and water tested prior to shipment.


### Available Options

- GC** Powder coated finish on galvanized Pipe and Fittings:. Available colors include orange, yellow, red and green.
- FC20** Regulates shower flow rate to 20 GPM.
- GRN** Green ABS plastic shower head.
- YEL** Yellow ABS plastic shower head.
- SSH** Stainless steel shower head.
- AP275-100 Electric Light and Alarm Horn**  
Light and alarm horn summon assistance when eyewash or shower is activated.
- AP250-015 Modesty Curtain**  
Modesty curtain for mounting on safety station.
- TMV** AP3800 thermostatic mixing valve precisely blends hot and cold water to deliver warm (tepid) water as required by ANSI Z358.1-2014.



**SS902HS** Safety Station with Eyewash, Auxiliary Drench Hose, Stainless Steel Bowl

**SS902P-HS** Safety Station with Eyewash, Auxiliary Drench Hose, Plastic Bowl

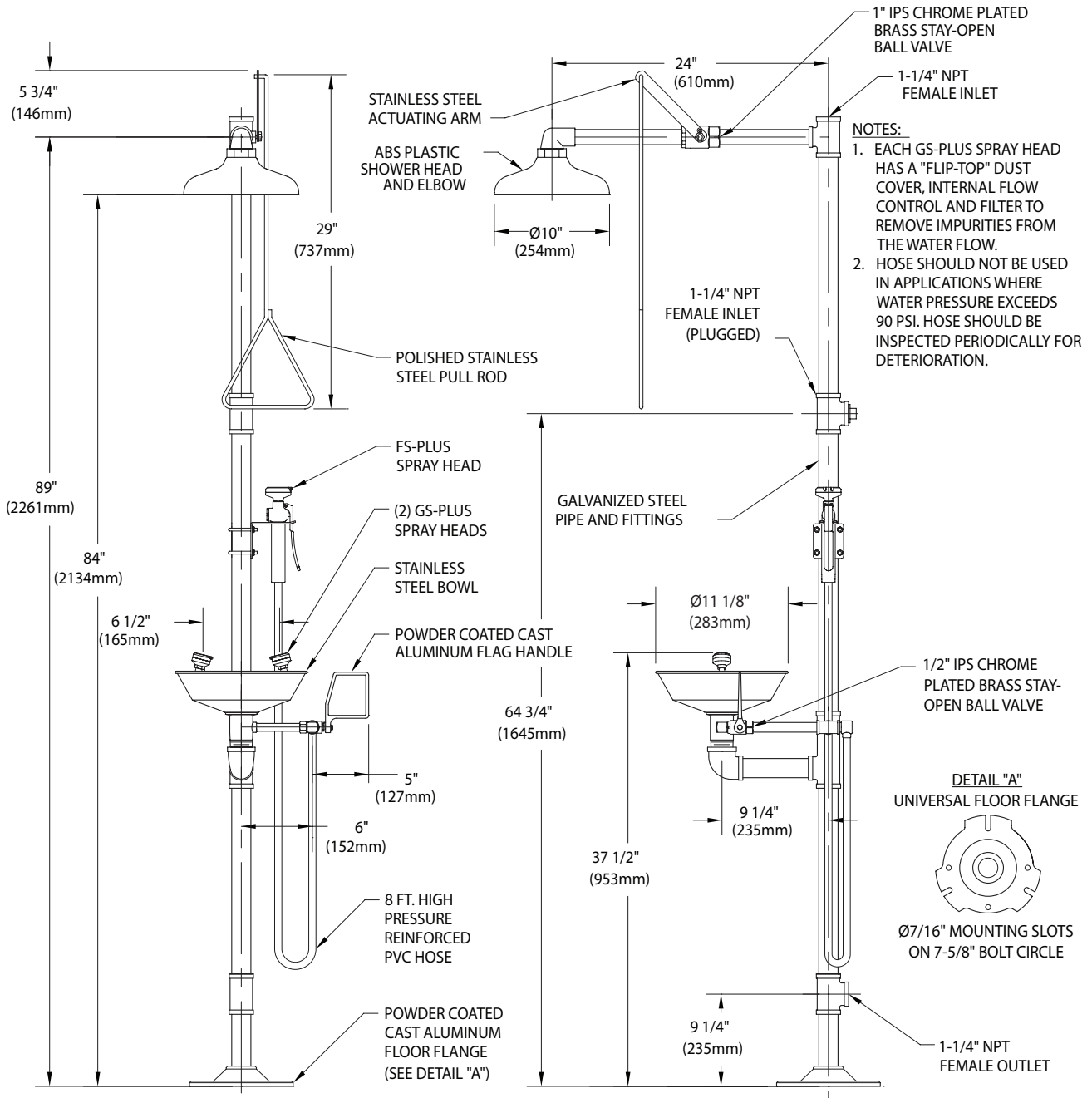

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are ± 1/4" (6mm). rev. 0615

Sign Included



WaterSaver  
701 W Erie St  
Chicago, IL 60654

312 666 5500 TELEPHONE  
312 666 5501 FACSIMILE  
wsflab.com

Listed 8116. Units have been tested to and comply with ANSI Z358.1-2014 and the Uniform Plumbing Code.



ETL Listed 101496. Units have been tested to and comply with ANSI Z358.1-2014.

