

**BT1614WSA** Mixer Tap, Bench Mounted



**Application:** Laboratory water mixer tap for hot and cold water.

**Mounting:** Bench or countertop mounted.

**Body:** Cast brass.

**Headworks:** WaterSaver self-contained compression valve cartridges with replaceable stainless steel seats.

**Handles:** Molded nylon hooded style handles with index discs. Handles and indexes are color coded and marked per EN 13792.

**Gooseneck:** 150mm (6") convertible rigid/swing gooseneck.

**Inlets:** (2) G1/2 male inlets on 180mm center dimension. Furnished with G1/2 mounting shanks, locknuts and washers.

**Outlet:** G3/8 female outlet with removable anti-splash serrated hose end.

**Quality Assurance:** Fitting is fully assembled and factory tested prior to shipment. 10 bar (150 PSI) test pressure and maximum working pressure.

**Ordering Information**

**Options**

<b>Valve Construction</b>	Compression valve cartridges	<input type="checkbox"/> BT1614WSA
	Ceramic disc valve cartridges	<input type="checkbox"/> BTC1614WSA
<b>Gooseneck</b>	200mm (8") gooseneck in place of 150mm (6")	<input type="checkbox"/> BT1614-8WSA
<b>Outlet Fitting</b>	BO055 aerator in place of hose end	<input type="checkbox"/> BT1614-55WSA
	BT118 pipe interrupter fitting with hose end	<input type="checkbox"/> BT1614-118WSA

**Product Finish**


<input type="checkbox"/> WHT	Gloss white powder coated finish	<input type="checkbox"/> PCL	Polished chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRY	Gloss light gray powder coated finish	<input type="checkbox"/> SCC	Satin chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRD	Gloss dark gray powder coated finish	<input type="checkbox"/> SNC	Satin nickel plated finish with clear epoxy coating
<input type="checkbox"/> STM	Starburst metallic powder coated finish		

**Note:** Refer to "Product Finishes" catalog for complete information on finishes.

**Media**

\_\_\_\_\_

BT1614WSA Mixer Tap, Bench Mounted

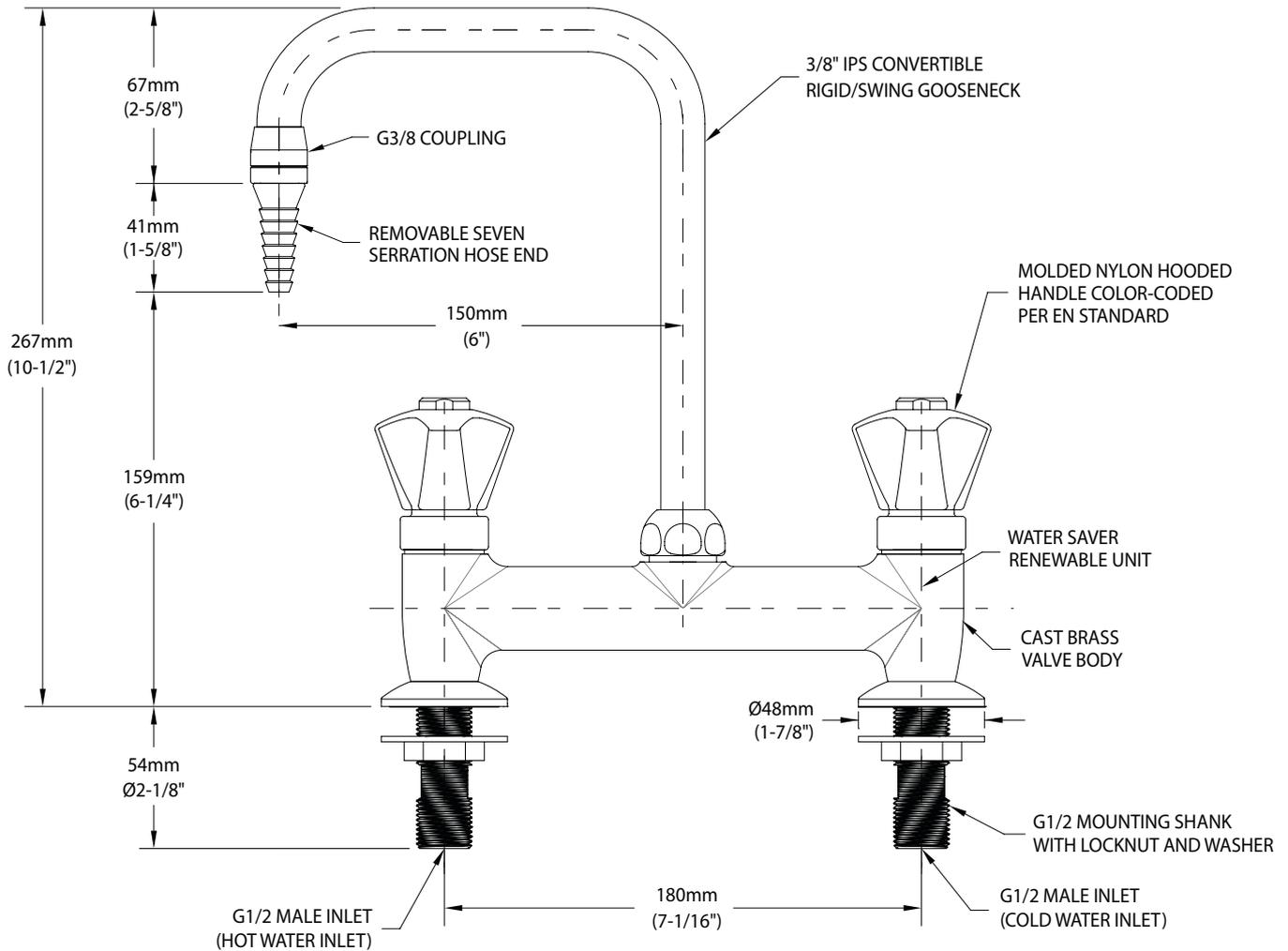

- NOTES:**
1. FITTING IS FULLY ASSEMBLED AND FACTORY TESTED PRIOR TO SHIPMENT.
  2. MAXIMUM WORKING PRESSURE IS 10 BAR (150 PSI).

Rev. 052318

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are  $\pm 6\text{mm}$  (1/4"). Rev. 071618