

- BT381** Double Water Tap, Bench Mounted, Two Valves at 180°
- BT382** Double Water Tap, Bench Mounted, Two Valves Side-by-Side



**Application:** Laboratory water tap with two outlets for cold water.

**Mounting:** Bench or countertop mounted with locator pins. On BT381, valves are mounted at 180° to each other. On BT382, valves are mounted side-by-side.

**Bodies:** Forged brass.

**Headworks:** WaterSaver self-contained compression valve cartridges with replaceable stainless steel seats.


**Handles:** Molded nylon hooded style handles with index discs. Handles and indexes are color coded and marked per EN 13792.

**Riser:** Brass riser with 150mm (6") spread to each valve.


**Inlet:** Furnished with G1/2 mounting shank, locknut and washer. G1/2 male inlet.

**Outlets:** G3/8 female outlets with removable anti-splash serrated hose ends.

**Quality Assurance:** Fitting is fully assembled and factory tested prior to shipment. 10 bar (150 PSI) test pressure and maximum working pressure.

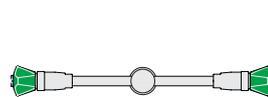

BT381

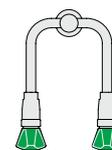

BT382

### Ordering Information

Options		BT381	BT382
<b>Valve Construction</b>	Compression valve cartridge	<input type="checkbox"/> BT381	<input type="checkbox"/> BT382
	Ceramic disc valve cartridge	<input type="checkbox"/> BTC381	<input type="checkbox"/> BTC382
<b>Riser</b>	Riser with 200mm (8") spread in place of 150mm (6")	<input type="checkbox"/> BT381-8	<input type="checkbox"/> BT382-8
<b>Outlet Fitting</b>	BO055 aerators in place of hose ends	<input type="checkbox"/> BT381-55	<input type="checkbox"/> BT382-55
	BT118 pipe interrupter fittings with hose ends	<input type="checkbox"/> BT381-118	<input type="checkbox"/> BT382-118

### Product Finish





<input type="checkbox"/> WHT	Gloss white powder coated finish	<input type="checkbox"/> PCL	Polished chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRY	Gloss light gray powder coated finish	<input type="checkbox"/> SCC	Satin chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRD	Gloss dark gray powder coated finish	<input type="checkbox"/> SNC	Satin nickel plated finish with clear epoxy coating
<input type="checkbox"/> STM	Starburst metallic powder coated finish		

**Note:** Refer to "Product Finishes" catalog for complete information on finishes.

### Media

\_\_\_\_\_

- BT381** Double Water Tap, Bench Mounted, Two Valves at 180°
- BT382** Double Water Tap, Bench Mounted, Two Valves Side-by-Side

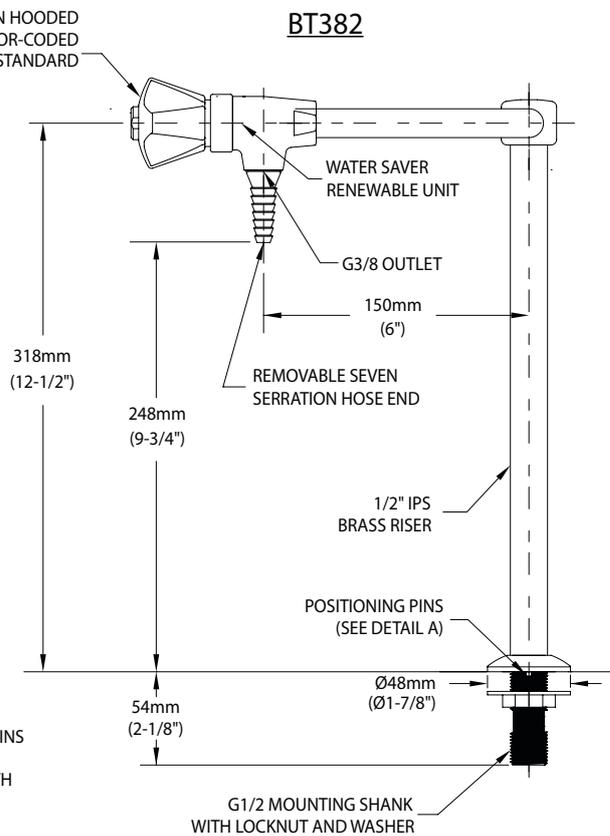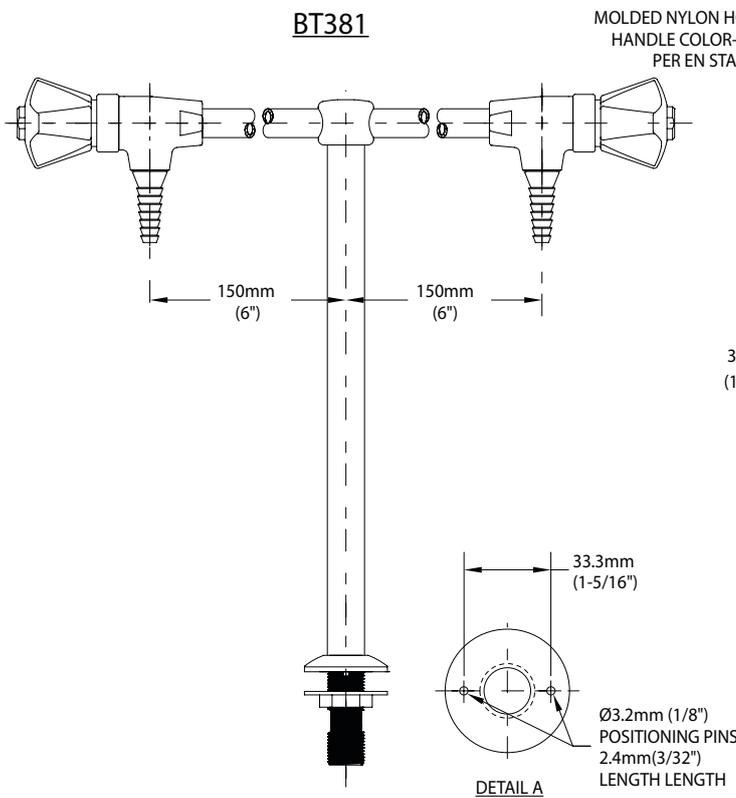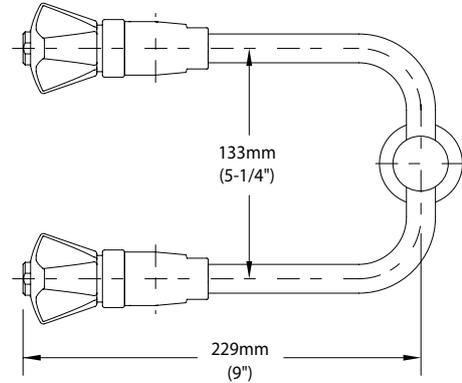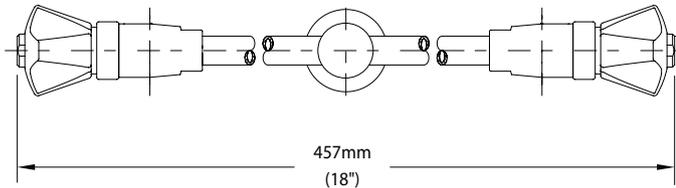

- NOTES:**
1. FITTING IS FULLY ASSEMBLED AND FACTORY TESTED PRIOR TO SHIPMENT.
  2. MAXIMUM WORKING PRESSURE IS 10 BAR (150 PSI).

Rev. 052318

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are ± 6mm (1/4"). Rev. 052118