

**BT385RS** Three-Way Water Tap, Bench Mounted



**Application:** Laboratory water tap with three outlets for cold water.

**Mounting:** Bench or countertop mounted with locator pins.

**Bodies:** Forged brass.

**Headworks:** WaterSaver self-contained compression valve cartridges with replaceable stainless steel seats.

**Handles:** Molded nylon hooded style handles with index discs. Handles and indexes are color coded and marked per EN 13792.

**Riser:** Brass riser with 150mm (6") spread to each side valve.

**Spout:** 150mm (6") convertible rigid/swing top spout.

**Inlet:** Furnished with G1/2 mounting shank, locknut and washer. G1/2 male inlet.

**Outlets:** G3/8 female outlets with removable anti-splash serrated hose ends.

**Quality Assurance:** Fitting is fully assembled and factory tested prior to shipment. 10 bar (150 PSI) test pressure and maximum working pressure.

**Ordering Information**

**Options**

<b>Valve Construction</b>	Compression valve cartridge	<input type="checkbox"/>	BT385RS
	Ceramic disc valve cartridge	<input type="checkbox"/>	BTC385RS
<b>Top Spout</b>	200mm (8") spout in place of 150mm (6")	<input type="checkbox"/>	BT385RS8
<b>Outlet Fitting</b>	BO055 aerator in place of hose end	<input type="checkbox"/>	BT385RS55
	BT118 pipe interrupter fitting with hose ends	<input type="checkbox"/>	BT385RS118

**Product Finish**


<input type="checkbox"/>	WHT	Gloss white powder coated finish	<input type="checkbox"/>	PCL	Polished chrome plated finish with clear epoxy coating
<input type="checkbox"/>	GRY	Gloss light gray powder coated finish	<input type="checkbox"/>	SCC	Satin chrome plated finish with clear epoxy coating
<input type="checkbox"/>	GRD	Gloss dark gray powder coated finish	<input type="checkbox"/>	SNC	Satin nickel plated finish with clear epoxy coating
<input type="checkbox"/>	STM	Starburst metallic powder coated finish			

**Note:** Refer to "Product Finishes" catalog for complete information on finishes.

**Media**

\_\_\_\_\_

BT385RS Three-Way Water Tap, Bench Mounted

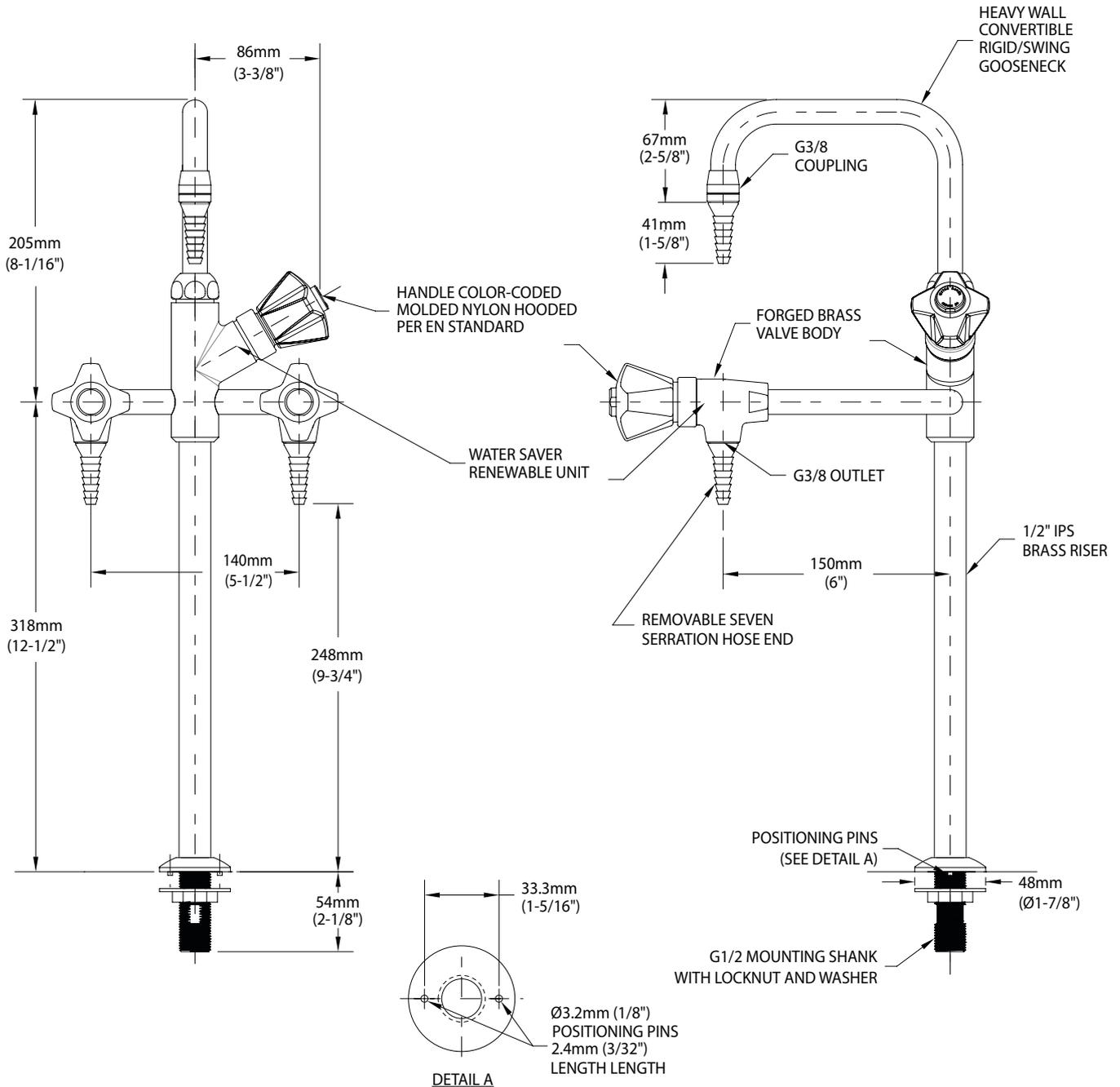

- NOTES:**
1. FITTING IS FULLY ASSEMBLED AND FACTORY TESTED PRIOR TO SHIPMENT.
  2. MAXIMUM WORKING PRESSURE IS 10 BAR (150 PSI).

Rev. 052318

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are ± 6mm (1/4"). Rev. 072618