

- BT681** Pure Water Tap, Bench Mounted, Self-Closing Control
- BT681LE** Pure Water Tap, Bench Mounted, Self-Closing Control, Lever Handle



Application: Laboratory tap for distilled, deionised or reverse osmosis water. All components in contact with pure water have an interior lining of pure tin.

Mounting: Bench or countertop mounted with locator pins.

Body: Forged brass.

Headwork: WaterSaver self-contained self-closing valve cartridge with replaceable stainless steel seat. Valve closes when handle is released.


Handle: Forged brass four-arm style handle (BT681) or forged brass lever handle (BT681LE), as ordered.

Gooseneck: 150mm (6") rigid gooseneck. Gooseneck is set in right hand configuration unless left hand or 180° is specified.

Inlet: Furnished with G1/2 mounting shank, locknut and washer. G1/2 male inlet.

Outlet: G3/8 female outlet with removable anti-splash serrated hose end.

Quality Assurance: Fitting is fully assembled and factory tested prior to shipment. 10 bar (150 PSI) test pressure and 7 bar (100 PSI) maximum working pressure.

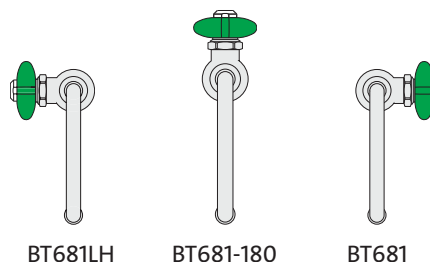

BT681LH BT681-180 BT681


Ordering Information			
Options			
	BT681	BT681LE	
Gooseneck	150mm (6") gooseneck	<input type="checkbox"/> BT681	<input type="checkbox"/> BT681LE
	200mm (8") gooseneck in place of 150mm (6")	<input type="checkbox"/> BT681-8	<input type="checkbox"/> BT681-8LE
Outlet Fitting	BO055 aerator in place of hose end	<input type="checkbox"/> BT681-55	<input type="checkbox"/> BT681-55LE
	BT118 pipe interrupter fitting with hose end	<input type="checkbox"/> BT681-118	<input type="checkbox"/> BT681-118LE

Product Finish			
<input type="checkbox"/> WHT	Gloss white powder coated finish	<input type="checkbox"/> PCL	Polished chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRY	Gloss light gray powder coated finish	<input type="checkbox"/> SCC	Satin chrome plated finish with clear epoxy coating
<input type="checkbox"/> GRD	Gloss dark gray powder coated finish	<input type="checkbox"/> SNC	Satin nickel plated finish with clear epoxy coating
<input type="checkbox"/> STM	Starburst metallic powder coated finish		

Note: Refer to "Product Finishes" catalog for complete information on finishes.

Media
<input type="checkbox"/> _____

- BT681** Pure Water Tap, Bench Mounted, Self-Closing Control
- BT681LE** Pure Water Tap, Bench Mounted, Self-Closing Control, Lever Handle

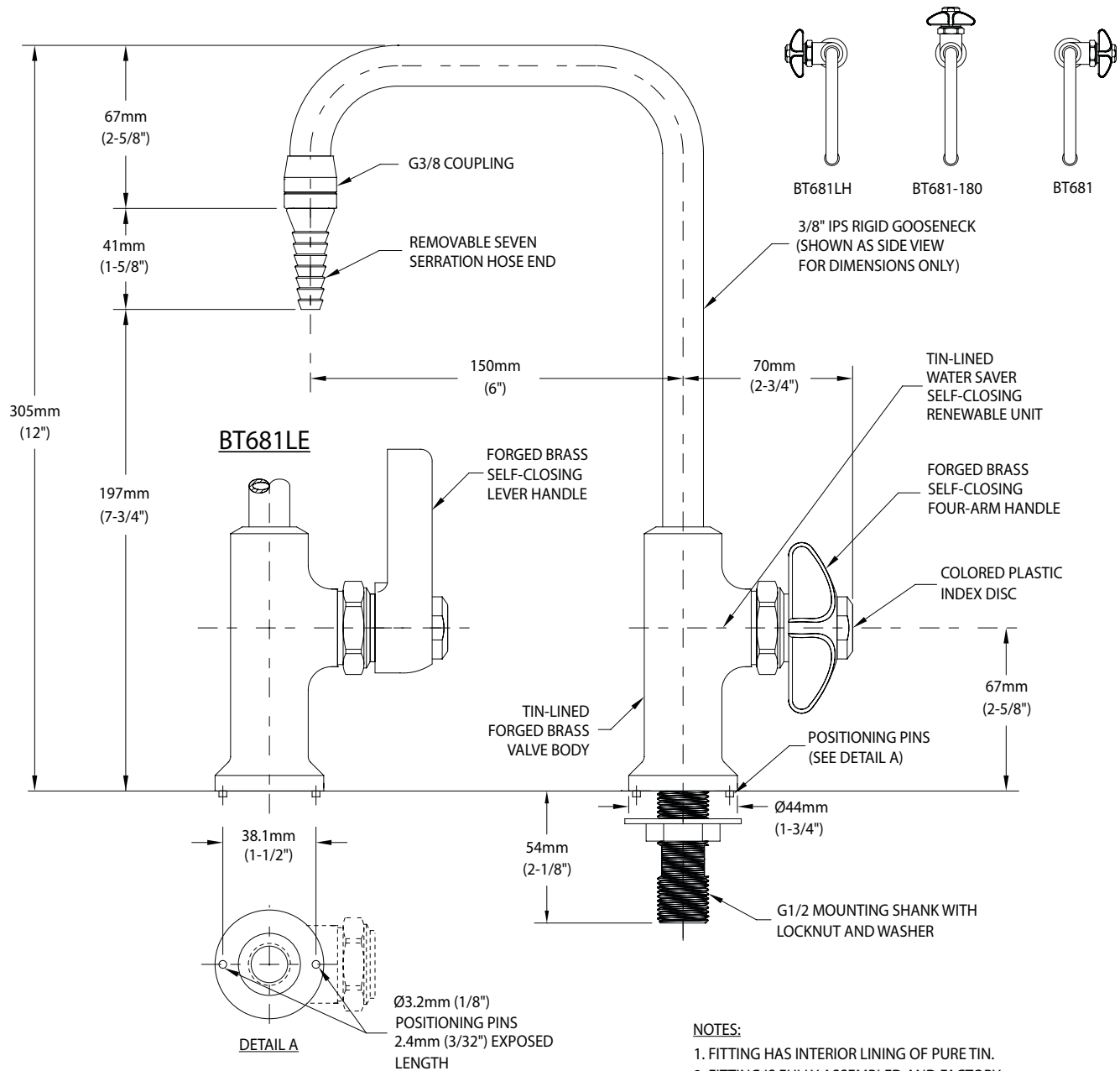

Rev. 031819

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are ± 6mm (1/4"). Rev. 031819