

BT691 Pure Water Tap, Bench Mounted, Manual Control

BT691BH Pure Water Tap, Bench Mounted, Manual Control, Wrist Blade Handle



Application: Laboratory tap for distilled, deionised or reverse osmosis water. All components in contact with pure water have interior lining of pure tin.

Mounting: Bench or countertop mounted with locator pins.

Body: Forged brass.

Headwork: WaterSaver self-contained compression valve cartridge with replaceable stainless steel seat.


Handle: Molded nylon hooded style handle (BT691) or forged brass wrist blade handle (BT691BH).

Gooseneck: 150mm (6") rigid gooseneck. Gooseneck is set in right hand configuration unless left hand or 180° is specified.

Inlet: Furnished with G1/2 mounting shank, locknut and washer. G1/2 male inlet.

Outlet: G3/8 female outlet with removable anti-splash serrated hose end.

Quality Assurance: Fitting is fully assembled and factory tested prior to shipment. 10 bar (150 PSI) test pressure and maximum working pressure.

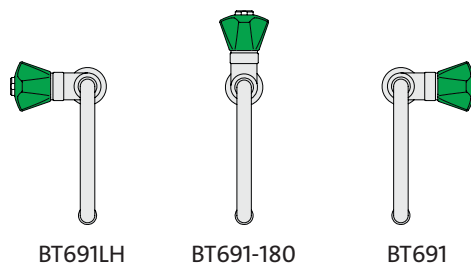

BT691LH

BT691-180

BT691

Ordering Information

| Options | | BT691 | BT691BH |
|-----------------------|--|------------------------------------|--------------------------------------|
| Gooseneck | 150mm (6") gooseneck | <input type="checkbox"/> BT691 | <input type="checkbox"/> BT691BH |
| | 200mm (8") gooseneck in place of 150mm (6") | <input type="checkbox"/> BT691-8 | <input type="checkbox"/> BT691-8BH |
| Outlet Fitting | BO055 aerator in place of hose end | <input type="checkbox"/> BT691-55 | <input type="checkbox"/> BT691-55BH |
| | BT118 pipe interrupter fitting with hose end | <input type="checkbox"/> BT691-118 | <input type="checkbox"/> BT691-118BH |

Product Finish

| | | | |
|------------------------------|---|------------------------------|--|
| <input type="checkbox"/> WHT | Gloss white powder coated finish | <input type="checkbox"/> PCL | Polished chrome plated finish with clear epoxy coating |
| <input type="checkbox"/> GRY | Gloss light gray powder coated finish | <input type="checkbox"/> SCC | Satin chrome plated finish with clear epoxy coating |
| <input type="checkbox"/> GRD | Gloss dark gray powder coated finish | <input type="checkbox"/> SNC | Satin nickel plated finish with clear epoxy coating |
| <input type="checkbox"/> STM | Starburst metallic powder coated finish | | |

Note: Refer to "Product Finishes" catalog for complete information on finishes.

Media
