
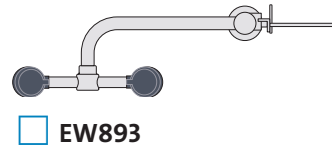


- ☐ **EW893** Eyewash, Deck Mounted 90° Swivel, All-Stainless Steel, Right Hand Mounting
- ☐ **EW893LH** Eyewash, Deck Mounted 90° Swivel, All-Stainless Steel, Left Hand Mounting

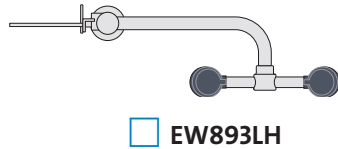

Application: All-stainless steel eyewash for mounting on counter next to sink. Unit is constructed entirely in stainless steel. Ideal for use in highly corrosive environments, cleanrooms and pharmaceutical applications. Spray heads swivel 90 degrees from storage to operational position. Water flow is activated by flag handle.

Spray Head Assembly: Two GS-Plus™ spray heads. Each head has a “flip top” dust cover, internal flow control and filter to remove impurities from the water flow.

Valve: 1/2" IPS stainless steel stay-open ball valve. Valve has stainless steel ball and PTFE seals.

Mounting: Unit may be mounted on right side (EW893) or left side (EW893LH) of sink. Furnished with locknut and washer for securing unit to countertop.

Construction: Type 316 stainless steel.

Supply: 1/2" NPT male inlet.


Sign: ANSI-compliant identification sign.

Quality Assurance: Unit is completely assembled and water tested prior to shipment.

Available Options

- ☐ **AP3600** Thermostatic mixing valve precisely blends hot and cold water to deliver warm (tepid) water as required by ANSI Z358.1-2014. Refer to “Tempering Valves” section for complete technical and product selection information.

- ☐ **EW893** Eyewash, Deck Mounted, 90° Swivel, All-Stainless Steel, Right Hand Mounting
- ☐ **EW893LH** Eyewash, Deck Mounted, 90° Swivel, All-Stainless Steel, Left Hand Mounting

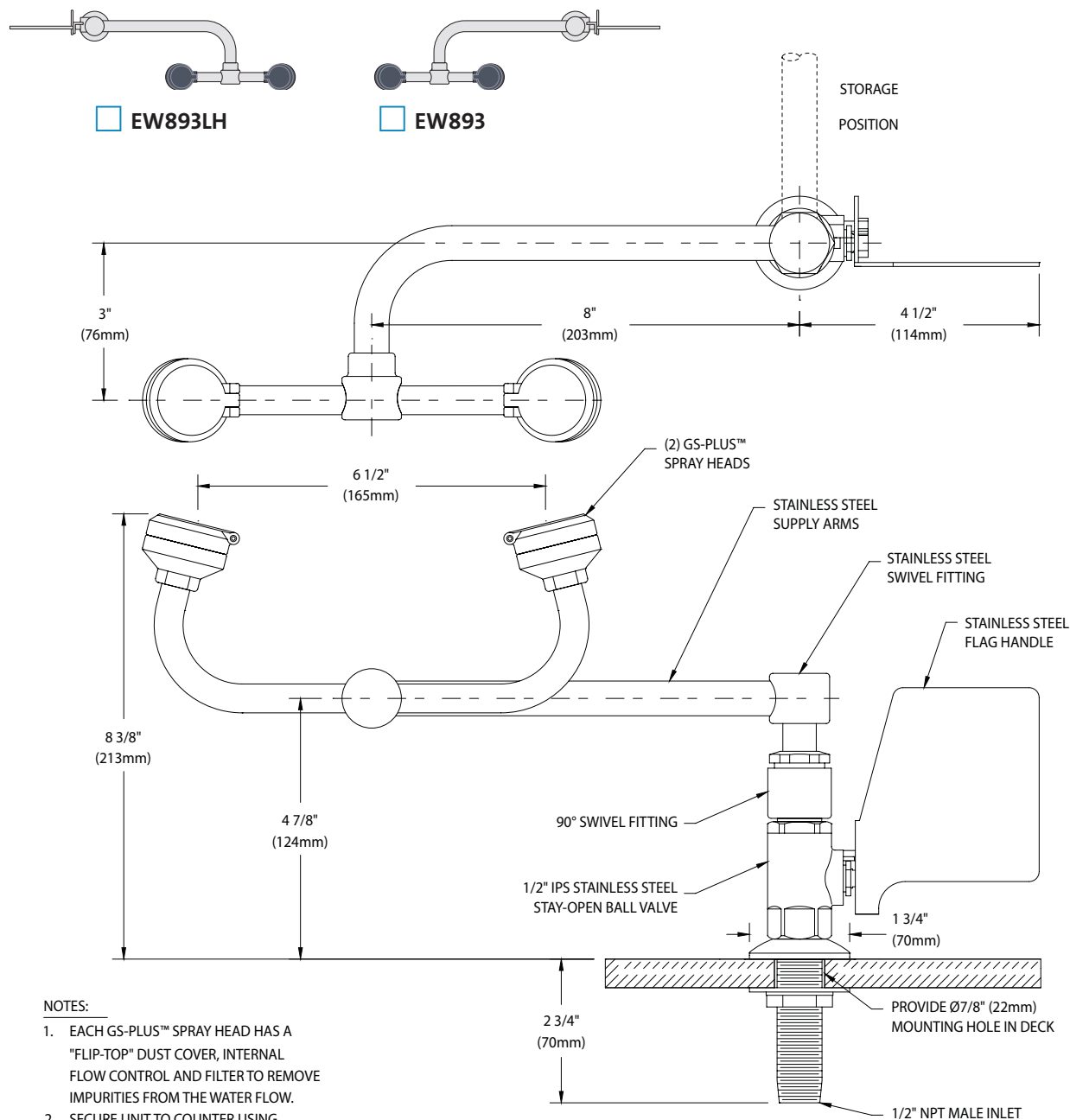

NOTES:

1. EACH GS-PLUS™ SPRAY HEAD HAS A "FLIP-TOP" DUST COVER, INTERNAL FLOW CONTROL AND FILTER TO REMOVE IMPURITIES FROM THE WATER FLOW.
2. SECURE UNIT TO COUNTER USING MOUNTING SHANK AND LOCKNUT.

THIS SPACE FOR ARCHITECT/ENGINEER APPROVAL

Sign Included



Due to continuing product improvement, the information contained in this document is subject to change without notice. All dimensions are $\pm 1/4"$ (6mm). rev. 011420

WaterSaver
701 W Erie St
Chicago, IL 60654

312 666 5500 TELEPHONE
312 666 5501 FACSIMILE
wsflab.com

Listed 8116. Units have been tested to and comply with ANSI Z358.1-2014 and the Uniform Plumbing Code.

